

Thermo Scientific Dionex ICS-4000
Integrated Capillary HPIC System



Always ready, better results
capillary high-pressure ion chromatography

Thermo
SCIENTIFIC

Integrated capillary ion chromatography system

The Dionex ICS-4000 HPIC System The World's First Dedicated Capillary High-Pressure Reagent-Free Ion Chromatography System

The Thermo Scientific™ Dionex™ ICS-4000 HPIC™ system combines a new level of resolution and speed with capillary ion chromatography, delivering high-pressure IC on demand with a system that is always on and always ready for the next analysis.

Benefits of Capillary HPIC

- Fast separations without compromising resolution
- True on-request IC analysis for 24/7 uptime with no equilibration
- HPIC performance for high chromatographic resolution
- Low eluent consumption and waste generation resulting in low cost of operation

Versatility

- Reagent-Free™ (RFIC™) system operation for isocratic and gradient separations with consistent results—day-to-day, lab-to-lab and operator-to operator
- Integrated system with greater performance and adaptability than a modular system
- Conductivity, electrochemical, and charge detection capabilities
- High performance system in a small footprint—occupies minimal bench space

Did You Know?

RFIC systems use deionized water to electrolytically generate precise eluent concentrations, and to regenerate suppressors, trap columns, and carbonate removal devices. This improves reproducibility, reduces labor, and improves detection limits.



Touchpad Display

- Displays key operating conditions
- Provides system set-up and control at the instrument

Dual-Stepper Motor Pump

- Precise and accurate capillary flow performance
- High-pressure, to support 4 μm particle columns and HPIC separations

Seal Wash (Option)

- Extra measure of reliability for pump seals where manually prepared eluents are used

Thermally Controlled Chromatography and Detection Zone

- Circulating air for cooling and heating
- Stable environment for column, suppressor, and detectors

Vacuum Degas

- In-line vacuum degassing of eluent
- Prevents eluent contamination from dissolved gases

Versatile Detector Options

- Conductivity or Electrochemical Detector
- New QD Charge Detector can also be configured with Conductivity Detector for complimentary detection

Tilt-Up Lid

- RFIC consumables located out-of-sight yet readily accessible
- Display tilts forward to allow local control while servicing



Reagent-Free Ion Chromatography

- Eluents generated electrolytically through software-control
- Reagent use of only 15 mL of water per day with continuous operation
- Eluent generator cartridge lasts up to 18 months of continuous operation under typical conditions

Eluent Shut-Off Valve

- Electrically controlled valve provides positive shut-off of source for easy pump servicing

Innovative Capillary IC Cube Design

- Dionex IC Cube™ design provides convenient plumbing of columns, suppressors, and other flow components
- Includes precut tubing and color-coded connections for ease of configuration and plumbing
- Heating for the column cartridge

Space-Saving Footprint

- Integrated system delivers dedicated performance in a small footprint

Versatile detection capabilities

The Dionex ICS-4000 Capillary HPLC system offers versatile detection options, from electrochemical detection to suppressed conductivity detection to our newest innovation, the QD charge detector, which is designed for detecting all ions.



CD Conductivity Detection

- Innovative, self-contained conductivity detector
- Supported by industry leading self-regenerating suppressor technology
- Optimized for capillary flow rates

QD Charge Detection

- New detector used in combination with CD detection for peak purity and identification
- Simplified calibration with only one standard
- Provides an equivalent signal for all ions of the same charge and molar concentration
- High linear response for weakly and strongly disassociated ions such as amines, organic acids, and silicate


ED Electrochemical Detection

- Sensitive and selective detection of cyanide, azide, biogenic amines and carbohydrates
- Innovative palladium hydrogen (PdH) reference electrode for superior reproducibility and long lifetime




High resolution and speed


With over 20 different column chemistries now available for capillary IC, the Dionex ICS-4000 HPIC system delivers state-of-the-art chromatographic performance for a wide range of applications.

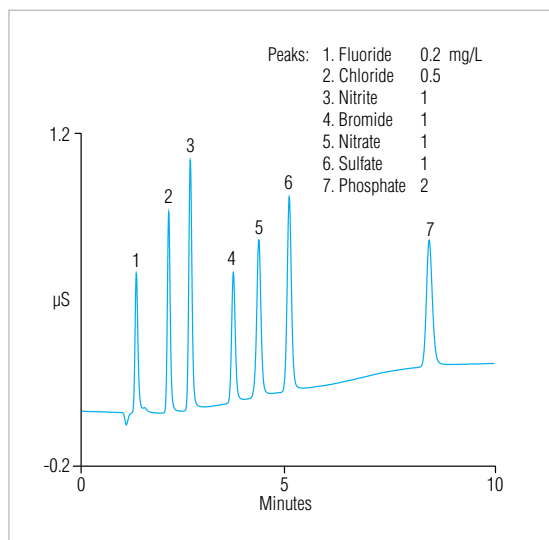


Gradient separation of inorganic anions in water using a Thermo Scientific Dionex IonPac™ AS18-4µm column and conductivity detection. The high-pressure capillary system and programmable electrolytic gradient lets you choose the proper balance between resolution and speed.

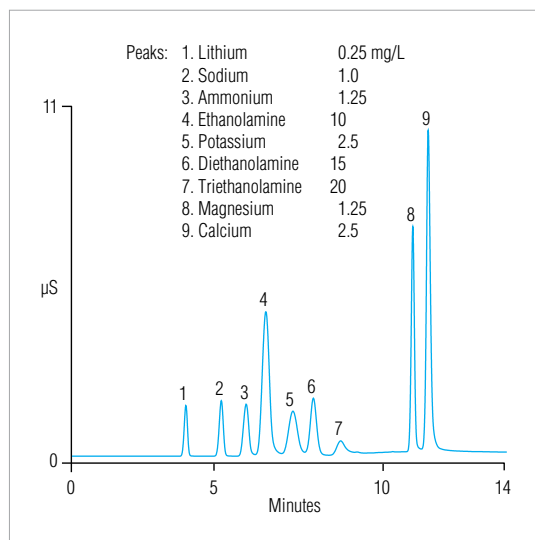

Gradient separation of ethanolamines in water using a Dionex IonPac CS19-4µm capillary column and conductivity detection. The high-pressure system and 4 µm column enables high resolution in half the time of a standard system.

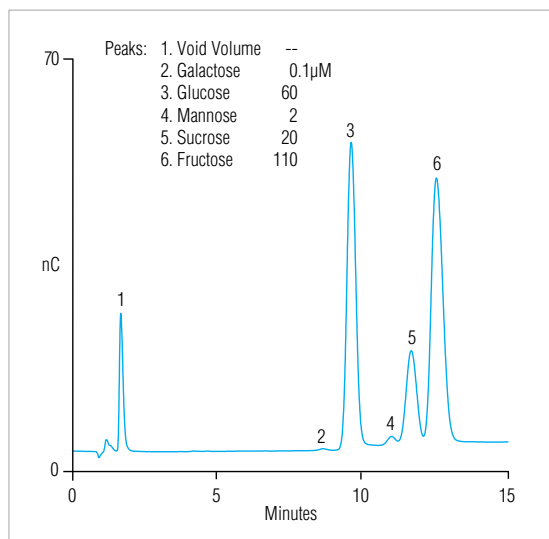

Determination of monosaccharides and disaccharides in apple cider using a Thermo Scientific Dionex CarboPac™ PA20 column and electrochemical detection. A wide range of concentrations can be separated and quantified.

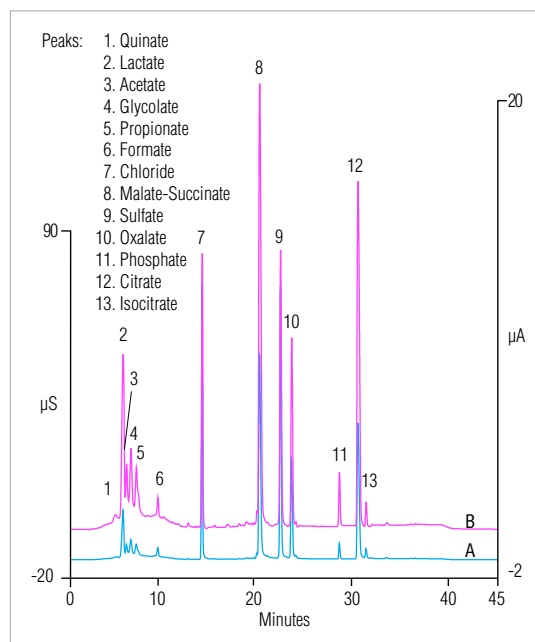

Comparison of A) suppressed conductivity and B) charge detection analysis of organic acids in guava juice, separated using a Dionex IonPac AS11-HC-4µm column and a KOH gradient.


Reliable sample handling

Versatile Sample Preparation and Automation Solutions

Dionex autosamplers integrate seamlessly with the Dionex ICS-4000 system, providing effortless automation and advanced sample preparation capabilities. All wetted surfaces are nonmetallic and chemically inert, protecting columns and sensitive samples from contamination.

Dionex AS-DV

- Affordable automation
- Precise mechanics for high reproducibility
- Sample filters integrated into vial caps
- Optional high-pressure valve for automated sample preparation
- Random access and sample overlap capabilities

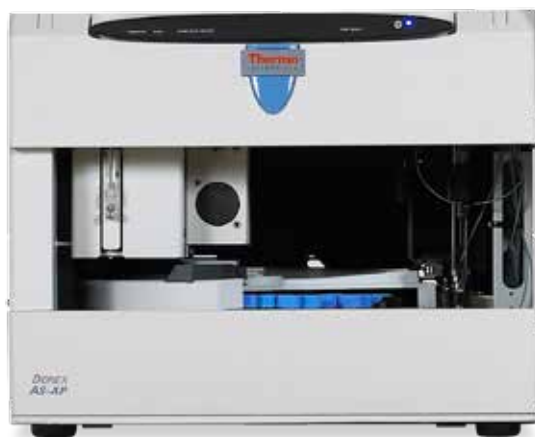




Dionex AS-AP

- Advanced performance autosampler with capillary mode and simultaneous and sequential injection support
- Automated sample preparation (dilution, concentration, matrix elimination, precolumn derivatization)
- Temperature control for thermally labile samples
- Fast injection cycle time, < 30 s per 5 μ L
- Flexible formats—10 mL vials to 0.45 mL well plates
- Optional Sample Conductivity and pH Accessory with conditionals


Simply intelligent software

Powerful Capabilities and Outstanding Usability

The natural complement to your IC system is Thermo Scientific™ Dionex™ Chromeleon™ software, the world's leading chromatography data system. With its inviting user interface and time-saving innovations, Chromeleon 7 software streamlines your path from samples to results. You'll find all the tools you need to tackle your toughest challenges without getting distracted from your goals. Get both rich functionality and outstanding usability with Chromeleon – the Simply Intelligent™ chromatography software.

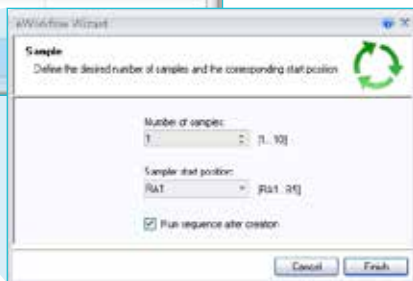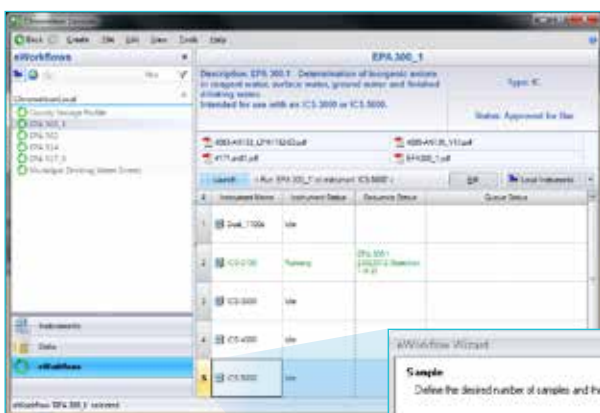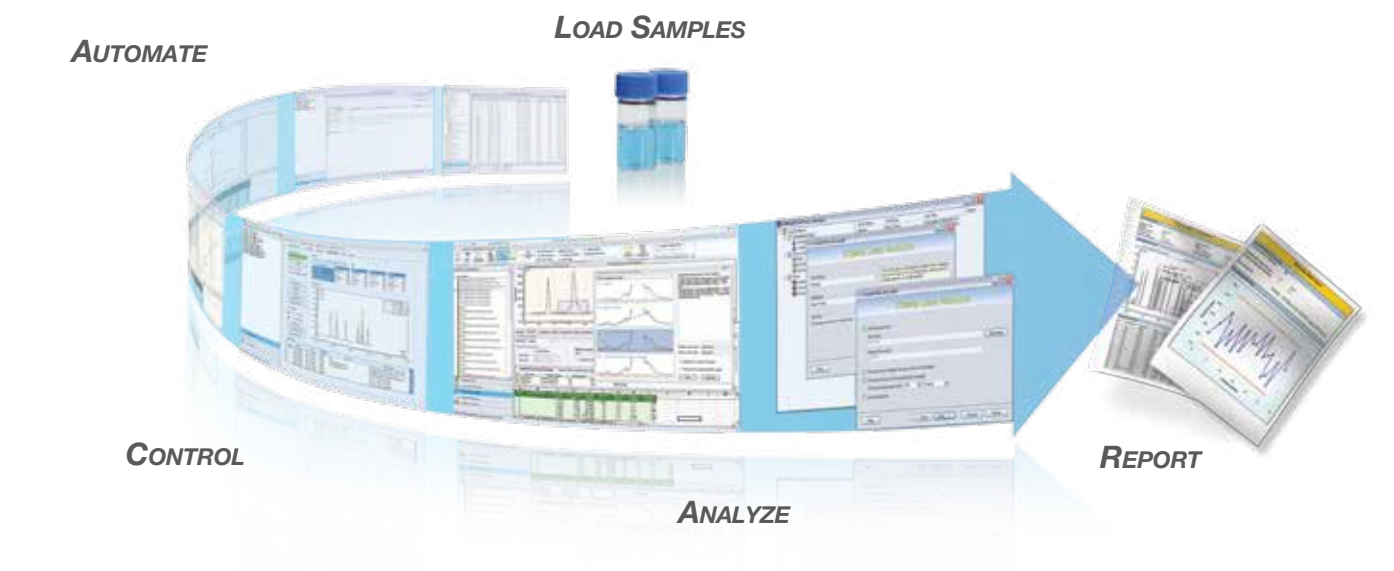

Gain an instant and intuitive understanding of your data at a glance with the Chromeleon 7 software's unique MiniPlots and Interactive Charts.

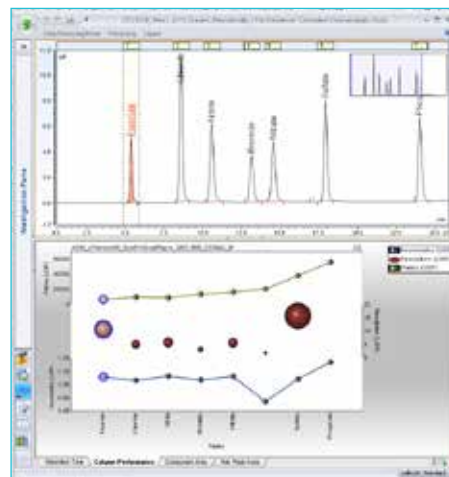

Discover the Dionex ICS-4000 system and HPIC

The Dionex ICS-4000 Capillary HPIC system is the next generation dedicated ion chromatography system, offering the ability to use 0.2–0.6 mm diameter column sets at high pressure. The high-pressure capability offers the flexibility to use small-particle columns for higher resolution or higher flow rates with standard or short columns for fast results.

The reduced eluent flow rate of capillary IC systems reduces eluent consumption and produces less waste. The Dionex IC Cube uses cartridge based consumables with color coded fittings, that are also compatible with alternate systems. Additionally, with an array of detectors and options, the ICS-4000 add up to a adds up to be a very versatile system.

For more information or to place an order, contact the Thermo Scientific Dionex Products office nearest you or your local distributor. Phone numbers and addresses for worldwide subsidiaries can be found in the About Us section of www.thermoscientific.com/dionex.

The Value of Support and Service

Customer care begins after delivery—when you purchase any Thermo Scientific instrument, module, or system, you get a partner, not just a piece of equipment. We provide access to a wide range of quality post-installation programs and assistance; including expert customer training, comprehensive service agreements, validation and qualification services, application support, and timely, professional service by skilled, locally based representatives. Our service representatives are distributed around the world, but close to you, ready to provide professional assistance, whether you're a new customer or a longtime partner. Your success equals our success.

www.thermoscientific.com/dionex

©2012 Thermo Fisher Scientific Inc. All rights reserved. ISO is a trademark of the International Standards Organization. All other trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.



Thermo Scientific Dionex products are designed, developed, and manufactured under an ISO 9001 Quality System.

Australia +61 3 9757 4486
Austria +43 1 616 51 25
Benelux +31 20 683 9768
+32 3 353 42 94
Brazil +55 11 3731 5140

China +852 2428 3282
Denmark +45 36 36 90 90
France +33 1 39 30 01 10
Germany +49 6126 991 0
India +91 22 2764 2735

Ireland +353 1 644 0064
Italy +39 02 51 62 1267
Japan +81 6 6885 1213
Korea +82 2 3420 8600
Singapore +65 6289 1190

Sweden +46 8 473 3380
Switzerland +41 62 205 9966
Taiwan +886 2 8751 6655
UK +44 1276 691722
USA and Canada +847 295 7500

Thermo
SCIENTIFIC
Part of Thermo Fisher Scientific